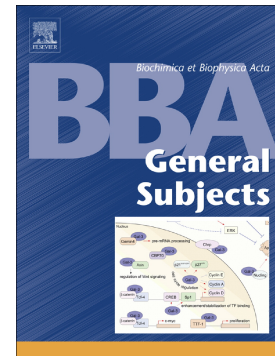


ReBaTSA: A simplified CeTSA protocol for studying recombinant mutant proteins in bacterial extracts

Maria Monticelli, Demi Marie Wright, Maria Vittoria Cubellis, Giuseppina Andreotti



PII: S0304-4165(23)00224-6

DOI: <https://doi.org/10.1016/j.bbagen.2023.130526>

Reference: BBAGEN 130526

To appear in: *BBA - General Subjects*

Received date: 9 October 2023

Revised date: 15 November 2023

Accepted date: 27 November 2023

Please cite this article as: M. Monticelli, D.M. Wright, M.V. Cubellis, et al., ReBaTSA: A simplified CeTSA protocol for studying recombinant mutant proteins in bacterial extracts, *BBA - General Subjects* (2023), <https://doi.org/10.1016/j.bbagen.2023.130526>

This is a PDF file of an article that has undergone enhancements after acceptance, such as the addition of a cover page and metadata, and formatting for readability, but it is not yet the definitive version of record. This version will undergo additional copyediting, typesetting and review before it is published in its final form, but we are providing this version to give early visibility of the article. Please note that, during the production process, errors may be discovered which could affect the content, and all legal disclaimers that apply to the journal pertain.

ReBaTSA: A SIMPLIFIED CETSA PROTOCOL FOR STUDYING RECOMBINANT MUTANT PROTEINS IN BACTERIAL EXTRACTS

Maria Monticelli^{1,2}, Demi Marie Wright^{1,3}, Maria Vittoria Cubellis^{1,2,4*} and Giuseppina Andreotti^{1*}

¹ Institute of Biomolecular Chemistry ICB, CNR, Via Campi Flegrei 34, 80078 Pozzuoli, Italy

² Department of Biology, University of Napoli “Federico II”, Complesso Universitario Monte Sant’Angelo, Via Cinthia, 80126 Napoli, Italy

³ Institute of Chemistry and Biochemistry, Freie Universität Berlin, D-14195 Berlin, Germany

⁴ Stazione Zoologica “Anton Dohrn”, Villa Comunale, Naples, Italy

* Correspondence: gandreotti@icb.cnr.it, cubellis@unina.it

KEYWORDS

CeTSA; ReBaTSA; Thermal Shift Assay; Protein stability; Phosphomannomutase-2; PMM2-CDG.

ABSTRACT

Introduction: The study of protein stability is crucial to biochemistry and relies on different methodologies. Recently, the Cellular Thermal Shift Assay has been introduced to study protein stability in whole cells.

Methods: We report a novel application of CeTSA named ReBaTSA. This Recombinant Bacterial TSA was performed using clear extracts from bacteria expressing a recombinant protein, incubated at different temperatures, centrifuged and analyzed via SDS-PAGE.

Results and Conclusions: We demonstrated the feasibility and reliability of this simplified approach. We validated the method using the protein phosphomannomutase-2 and its common mutants, which were compared in the presence or the absence of a known ligand.

MAIN TEXT

The melting temperature of a protein (T_M) is the temperature at which 50% of the molecules are in the native state. Its value depends on the difference in free energy between the native and denatured forms. The process must be reversible to correctly determine the melting temperature, the enthalpy of fusion, and other thermodynamic quantities. This condition is challenging because chemical reactions that lead to the irreversible denaturation of proteins accelerate as the temperature increases. Nevertheless, estimating T_M , even under non-reversible conditions can provide comparisons of numerous variants of the same protein.

Comparing the T_M of wild-type and variant proteins is advantageous for many reasons, ranging from biotechnological to medical applications. Most missense mutations associated with disease do not occur in the active site of proteins or other functionally important sites. Some mutations prevent folding; others cause a decrease in stability, manifesting in a decreased melting temperature¹. Evaluating the impact of a mutation on a protein is a complex but essential task in exomic analyses. Several *in silico* methods exist to assess whether the mutation is deleterious, and experimental datasets are used to validate these methods by measuring the melting temperature. Moreover,

analysing the effect of a mutation on a protein also allows for patient stratification and personalised therapies.

In some cases, unstable proteins exhibit a decreased intracellular concentration compared to their wt counterpart. Protein stability is dependent on several factors: environment, pH, ionic strength, and physical-chemical environment, as well as the presence of ligands. It has been demonstrated that in the presence of specific small molecules, corresponding mutated and unstable proteins exhibit an increased melting temperature. This apparent stabilisation leads to increase concentration – and consequently the total activity - of mutated proteins.

Cellular Thermal Shift Assay (CeTSA) is a versatile technique introduced in 2013 ^{2,3} as an alternative to Thermal Shift Assay (TSA ⁴). It measures the residual amount of protein after heating at a specific temperature. This value is sensitive to the environment, pH, ionic strength, protein-protein and protein-ligand interactions. CeTSA, as it is for TSA, can be used to plot melting profiles and estimate T_M , carrying out experiments at different temperatures and recording the ratio between the amount of protein after and before heating.

Isothermal CeTSA can also be used to assess the stabilizing effect of ligands. In this case, the temperature is kept constant and close to the T_M of the protein and the unfolded fraction is measured upon incubation with or without small molecules. Isothermal dose response fingerprint-CeTSA can be carried out by varying a small molecule's concentration (ITDRF-CeTSA) ^{5,6}.

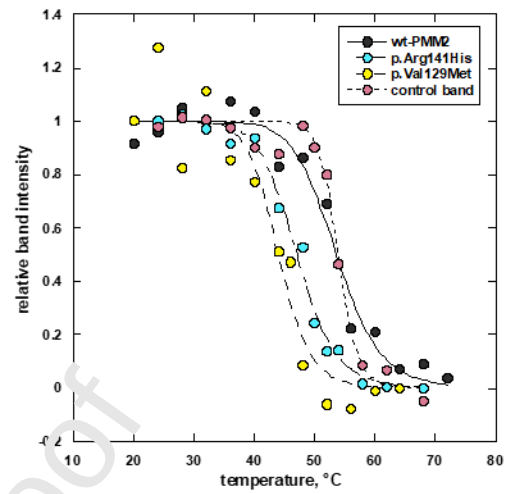
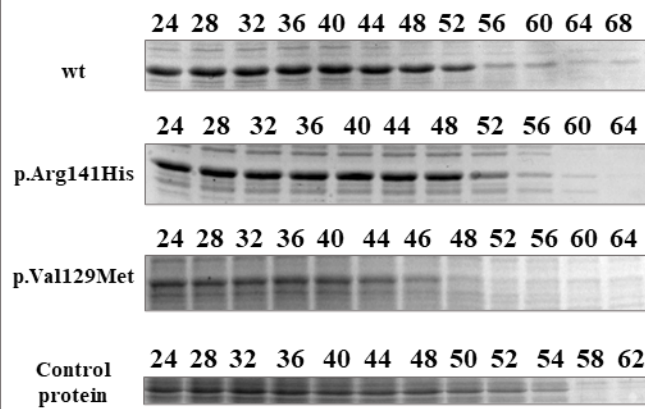
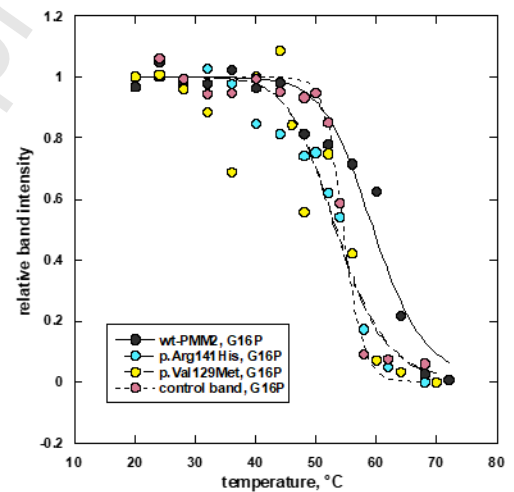
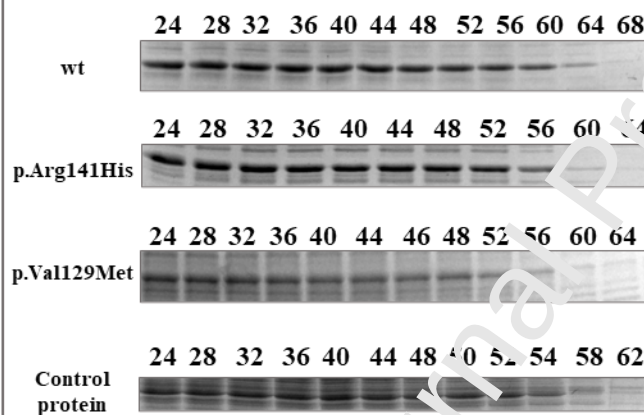
The novelty of CeTSA is related to its application to proteins in complex mixtures, such as cellular lysates or tissue homogenates. The original method required an immunoblot for the detection, and the latest versions used mass-spectrometry (MS CeTSA).

In this brief application note, we describe a simplified version of CeTSA on raw extracts from transformed bacteria that hyper-express the protein of interest, named ReBaTSA (Recombinant Bacterial TSA). It does not require protein purification, immunoblotting, or mass-spectrometry detection and can be adapted to isothermal analysis. The simplified method described below allows for comparison of wild-type protein with its unstable variants with or without ligands. The difficulty of purifying and preserving relatively unstable proteins should not be underestimated when evaluating the method's usefulness. Results obtained by ReBaTSA are comparable with those obtained by TSA. Our method can also be applied to isothermal CeTSA.

We tested ReBaTSA on the human enzyme phosphomannomutase-2 (PMM2). PMM2 converts mannose-6-phosphate (M6P) to mannose-1-phosphate (M1P) in the early glycosylation steps ⁷. It requires glucose-1,6-bisphosphate (G16P) for its activity ⁸. More than 115 mutations in phosphomannomutase-2 (PMM2) that cause congenital disorder of glycosylation (PMM2-CDG) have been described ⁹. Here, we focus our attention on two mutants of interest, namely p.Arg141His and p.Val129Met. p.Arg141His is the first most common pathogenic mutant worldwide ¹⁰. Arg141 occurs in the active site; hence, the variant is not active and has never been observed in homozygosis. It is often found in patients associated with hypomorphic mutants with which it forms heterodimers ¹¹. p.Val129Met is common in Italian patients and is often heterozygous with p.Arg141His ¹².

Wild-type and mutant PMM2 were expressed in *E. coli* as previously described ¹². Upon lysis, aliquots of raw protein extract were incubated at different temperatures in the 20-74°C range for 4 minutes each ¹³, then immediately cooled on ice. Samples were centrifuged (18400 x g) for 15 minutes at 4°, allowing the separation of pellets containing insoluble proteins and supernatants containing the soluble fraction. Supernatants were collected and analyzed via SDS-PAGE followed

by Coomassie staining (Fig. 1A). A quantitative analysis of the SDS-PAGE bands performed using Image Lab™ (Bio-Rad) allowed the determination of melting curves, where the band intensities of the target proteins were plotted versus the incubation temperatures. The T_M were then mathematically evaluated, applying the sigmoidal fit with KaleidaGraph™.

A**B****C**

PROTEIN	TSA*		ReBaTSA	
	$T_M, ^\circ\text{C}$		$T_M, ^\circ\text{C}$	
	-	+ G16P 0.5mM	-	+ G16P 0.5mM
wt-PMM2	53.3±0.8	59.9±3.8	51.8±1.7	57.9±1.5
p.Arg141His, PMM2	49.5±0.5	53.0±0.8	48.3±1.8	53.7±0.6
p.Val129Met, PMM2	47.7±0.6	58.1±1.3	45.2±1.6	52.7±1.0
Control protein	-	-	54.2±0.4	53.2±1.9

* Mean from data available in literature. Details in the caption.

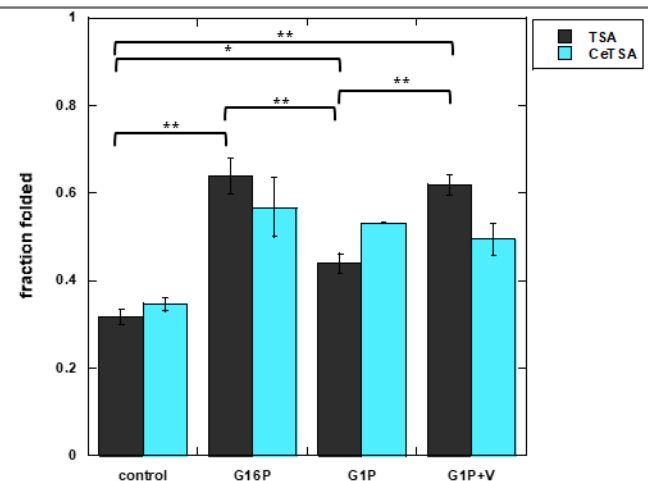
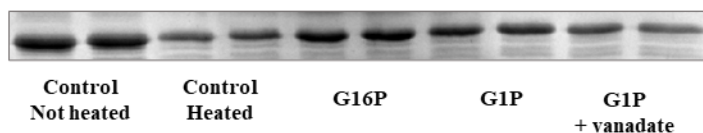
D

Figure 1. (A) Heat-induced melting profiles of *E. coli* extract expressing wt or mutant PMM2. The clear extracts in 50 mM Tris HCl pH 7.5, 1 mM MgCl₂, 0.1 mM DTT, were diluted up to 0.3 mg/mL and aliquots were incubated at different temperatures in the 20-74°C range for 4 minutes each, then immediately cooled on ice. Centrifuge at 18400 x g for 15 minutes at 4°C allowed pelletizing denatured proteins, and supernatants containing the folded fraction were collected and analyzed via SDS-PAGE. The Coomassie staining is showed on the left, quantitative analysis performed with the software ImageLab™ is shown on the right. (B) Heat-induced melting profiles of *E. coli* extract expressing wt or mutant PMM2 in the presence of 0.5 mM G16P. (C) The table summarizes the results (n=2). The mean of published results by fluorescence-based Thermal Shift Assay (TSA) - obtained with the purified wt or mutant PMM2 in the same experimental conditions - are shown for comparison^{7,8,11,12,14,15}. (D) Isothermal ReBaTSA. *E. coli* extract expressing wt-PMM2, was used and the effects of different ligands (published in [2], G16P, G1P, G1P + vanadate 0.5 mM each) were compared. The incubation was conducted at 54°C. Quantitative analysis of bands highlighted the different binding of the ligands considered (Tukey's HSD, * p < 0.05, ** p. < 0.01, n=2). For a comparison, data extracted from Fig.11 in Andreotti et al.¹⁵ - obtained by TSA using purified wt-PMM2 under the same experimental conditions - are shown.

Figure 1A shows the residual soluble target protein fraction compared to the untreated sample plotted versus the incubation temperature. T_M of wt-PMM2, p.Arg141His and p.Val129Met are 51.8±1.7°C, 48.3±1.8°C and 45.2±1.6°C respectively. In our experiments, the T_M perfectly aligned with those evaluated using fluorescence-based Thermal Shift Assay (TSA) on purified proteins^{7,8,11,12,14,15}.

The coincidence between the obtained results with ReBaTSA and TSA is also observed when the experiment is carried out in the presence of a ligand. Shift induced by G16P binding increased the T_M by 6°C for wt-PMM2, 5°C for p.Arg141His, and more than 7°C for p.Val129Met (Fig. 1 B).

The experiments were performed in analogous experimental conditions respect to published results, and G16P was used at a concentration 0.5 mM, as previously described. Figure 1C summarizes the results and shows the comparison with published data^{12,8,15,11,7,14}. Data from the papers above were extracted for comparison, and the mean was calculated.

ReBaTSA consents to use a non-specific abundant protein band as an internal control. The T_M of this protein did not change in the presence of G16P, thereby excluding a non-specific effect of the ligand (Fig. 1C).

Determining the melting temperature (T_M) of a protein is a preliminary step to carry out isothermal CeTSA, a simplified method to test the effect of ligands incubating the protein with or without the molecule of interest at a temperature close to the melting point.

Wt-PMM2 was incubated at a fixed temperature to compare various ligands. G16P, glucose-1-phosphate (G1P) and G1P plus vanadate were selected based on previous investigation¹⁵. The stabilizations were recorded at 54°C. The ratio between the residual amount of protein and its total amount was calculated by the bands' intensity and compared with the values obtained by melting curves previously carried out with the same protein and the same ligands by TSA (Fig 1D).

The IsoThermal Dose Response Fingerprint-CeTSA (ITDRF-CeTSA), dedicated to analyzing dose-response ligand binding, exhibits the development of this approach. We performed ITDRF-ReBaTSA on our target proteins using G16P. According to the melting profiles (Fig. 1A), we used a varying incubation temperature for each mutant: 53°C for p.Arg141His and p.Val129Met, 57°C for wt-PMM2. The results are shown in Supplementary Figure 1. As already mentioned, p.Arg141His is

an inactive variant of PMM2. For this reason, it was not possible to evaluate the affinity for G16P by enzymatic assays. Our experiments demonstrate that G16P binds to the inactive protein with an affinity comparable to that of wild-type protein, at least at high temperature.

In conclusion, this report demonstrated the feasibility and robustness of the ReBaTSA approach. The limitations associated with ReBaTSA protocol are mainly related to the recombinant protein expression. It is necessary that the target protein is expressed at quite high levels, to use Coomassie staining rather than immunoblot. Also, the effects mediated by relevant post-translational modifications would not be recorded, representing false negatives. Another restraint comes from the limited permeability of bacterial cells that hampers the possibility to test ligands using whole cells. Besides this, the study of protein complexes is not suitable with this approach. Despite its limitations, we believe it is a valuable biochemical tool for comparing wt and mutant variations of the target protein. This is an ever-expanding field since the study of mutated proteins covers a wide range of biotechnological applications, ranging from biomedical fields to bio-productive processes.

ABBREVIATIONS

CeTSA: Cellular Thermal Shift Assay

G16P: Glucose-1,6-Bisphosphate

G1P: Glucose-1-Phosphate

ITDRF-CeTSA: Isothermal dose response fingerprint-CeTSA

M1P: Mannose-1-Phosphate

M6P: Mannose-6-Phosphate

ReBaTSA: Recombinant Bacterial TSA

T_M : Melting temperature

TSA: Fluorescence-based Thermal Shift Assay

AUTHOR CONTRIBUTIONS

Conceptualization MVC and GA; Formal analysis MM, DMW, MVC and GA; Funding acquisition MVC; Investigation MM and DMW; Methodology GA; Resources MVC and GA; Software GA; Supervision MVC and GA; Validation MM, DMW and GA; Visualization MM, DMW and GA; Roles/Writing - original draft MM and GA; and Writing - review & editing DMW and MVC.

ACKNOWLEDGEMENTS

We acknowledge the National Research Council of Italy, Joint Bilateral Agreement CNR/Slovak Academy of Sciences, Biennial program 2023-24.

FUNDING

We acknowledge the funding by the Italian Ministry of University and Research PRIN 2022B2N2BY.

References

1. DePristo MA, Weinreich DM, Hartl DL. Missense meanderings in sequence space: A biophysical view of protein evolution. *Nat Rev Genet.* 2005;6(9). doi:10.1038/nrg1672
2. Jafari R, Almqvist H, Axelsson H, et al. The cellular thermal shift assay for evaluating drug target interactions in cells. *Nat Protoc.* 2014;9(9). doi:10.1038/nprot.2014.138
3. Tolvanen TA. Current Advances in CETSA. *Front Mol Biosci.* 2022;9:866764. doi:10.3389/fmolb.2022.866764
4. Andreotti G, Monticelli M, Cubellis MV. Looking for protein stabilizing drugs with thermal shift assay. *Drug Test Anal.* 2015;7(9). doi:10.1002/dta.1793
5. Ishii T, Okai T, Iwatani-Yoshihara M, et al. CETSA quantitatively verifies in vivo target engagement of novel RIPK1 inhibitors in various bi-specimens. *Sci Rep.* 2017;7(1):13000. doi:10.1038/s41598-017-12513-1
6. Friman T. Mass spectrometry-based Cellular Thermal Shift Assay (CETSA®) for target deconvolution in phenotypic drug discovery. *Bioorg Med Chem.* 2020;28(1):115174. doi:10.1016/j.bmc.2019.115174
7. Andreotti G, Pedone E, Giordano A, Cubellis MV. Biochemical phenotype of a common disease-causing mutation and a possible therapeutic approach for the phosphomannomutase 2-associated disorder of glycosylation. *Mol Genet Genomic Med.* 2013;1(1). doi:10.1002/mgg3.3
8. Citro V, Cimmaruta C, Liguori L, Viscido G, Cubellis MV, Andreotti G. A mutant of phosphomannomutase1 retains full enzymatic activity, but is not activated by IMP: Possible implications for the disease PMM2-CDG. *PLoS ONE.* 2017;12(12). doi:10.1371/journal.pone.0189629
9. PMM2 gene: MedlinePlus Genetics. Accessed October 4, 2023. <https://medlineplus.gov/genetics/gene/pmm2/>
10. Chen S, Francioli LC, Goodrich JK, et al. A genome-wide mutational constraint map quantified from variation in 76,156 human genomes. Published online October 10, 2022:2022.03.20.485034. doi:10.1101/2022.03.20.485034
11. Andreotti G, Monti MC, Citro V, Cubellis MV. Heterodimerization of two pathological mutants enhances the activity of human phosphomannomutase2. *PLoS ONE.* 2015;10(10). doi:10.1371/journal.pone.0139882
12. Citro V, Cimmaruta C, Monticelli M, et al. The analysis of variants in the general population reveals that PMM2 is extremely tolerant to missense mutations and that diagnosis of PMM2-CDG can benefit from the identification of modifiers. *Int J Mol Sci.* 2018;19(8). doi:10.3390/ijms19082218

13. Seashore-Ludlow B, Axelsson H, Lundbäck T. Perspective on CETSA Literature: Toward More Quantitative Data Interpretation. *SLAS Discov.* 2020;25(2):118-126. doi:10.1177/2472555219884524
14. Monticelli M, Liguori L, Allocca M, Andreotti G, Cubellis MV. β -Glucose-1,6-Bisphosphate stabilizes pathological phosphomannomutase2 mutants in vitro and represents a lead compound to develop pharmacological chaperones for the most common disorder of glycosylation, PMM2-CDG. *Int J Mol Sci.* 2019;20(17). doi:10.3390/ijms20174164
15. Andreotti G, De Vaca IC, Poziello A, Monti MC, Guallar V, Cubellis MV. Conformational response to ligand binding in Phosphomannomutase2: Insights into inborn glycosylation disorder. *J Biol Chem.* 2014;289(50). doi:10.1074/jbc.M114.586362

AUTHOR CONTRIBUTIONS

Conceptualization MVC and GA; Formal analysis MM, DMW, MVC and GA; Funding acquisition MVC; Investigation MM and DMW; Methodology GA; Resources MVC and GA; Software GA; Supervision MVC and GA; Validation MM, DMW and GA; Visualization MM, DMW and GA; Roles/Writing - original draft MM and GA; and Writing - review & editing DMW and MVC.

Declaration of interests

The authors declare that they have no known competing financial interests or personal relationships that could have appeared to influence the work reported in this paper.

HIGHLIGHTS

- Protein stability is essential to protein function, and missense destabilising mutations underlie many human diseases.
- The study of protein stability by biochemical and biophysical methodologies typically relies on protein expression and purification.
- Cellular Thermal Shift Assay (CeTSA) is an innovative technique for studying protein stability in whole cells. Since its first description, many applications have been described that require the identification of the protein of interest by immunoblots or mass-spectrometry.
- We show a proof-of-concept application of CeTSA to raw protein extracts for fast mutant comparison and patient stratification.
- This simplified approach bypasses protein purification and complex and expensive visualisation (i.e., immunoblot or mass spectrometry). The technique is specific and reliable.